

New!

JOBST® Anti-Embolism

**Designed to reduce
the risk of deep vein
thrombosis in the
recovering patient**

#1 physician recommended

JOBST®

Comfort, Health and Style!

JOBST® Anti-Embolism

Designed for the recovering patient

JOBST® Anti-Embolism Knee Highs and Thigh Highs are designed to reduce the risk of Deep Vein Thrombosis (DVT), commonly known as blood clots.

Specially designed as a Thrombo Embolic Deterrent stocking, JOBST Anti-Embolism applies gradient compression highest at the ankle and gradually decreases this pressure up the leg to improve blood flow. It is improved blood flow that helps reduce the risk for DVT.

JOBST Anti-Embolism Knee Highs and Thigh Highs are for patients who are resting in bed or otherwise immobile (non ambulatory).

JOBST Anti-Embolism Knee Highs and Thigh Highs are designed to be worn by both men and women.

Features & Benefits

Leg Therapy

- Accurate gradient compression therapy minimizes the risk of DVT (blood clots).

Comfortable

- Soft, breathable fabric is lightweight and comfortable.
- Thigh High style has a honeycomb silicone band to keep the product in place. This honeycomb design allows more air to circulate, promoting a cooler, more comfortable feeling.
- Latex-free formula helps avoid allergic reactions sometimes associated with other hosiery products.

Durable

- Reinforced heel is designed for long lasting wear.

Compliance

- The variety of sizes in knee and thigh high styles ensures a proper fit which increases patient compliance.

Styles

Knee High



Thigh High



Knee High and Thigh High are available for both Men and Women.

Color

White



JOBST® Anti-Embolism

Knee High and Thigh High... designed for all-day comfort!

Honeycomb silicone band keeps Thigh High in place

Fabric is lightweight and air-permeable

Reinforced heel for added durability and a better fit

All-day comfort knee band

Comfortable, closed toe, smooth, flat toe seam



Compare JOBST® Anti-Embolism to the competition:

Attribute	JOBST® Anti-Embolism	Kendal T.E.D.®	Activa® Anti-Embolism	Bell Horn®
Toe Seam Photo at same angle and distance				
Heel Photo at same angle and distance				
Thigh Band Photo at same angle and distance				
Knee Band Photo at same angle and distance				
Fabric Detailed photo of knit structure				

JOBST® Anti-Embolism


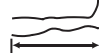
Measuring and Ordering Information:

KNEE HIGH

Size	Calf Circumference 	Length 		Catalog #	UPC Code
Small	9"-12"	Regular	14" - 17"	111470	0-35664-11470-1
		Long	17" - 20"	111471	0-35664-11471-8
Medium	12"-15"	Regular	14" - 17"	111472	0-35664-11472-5
		Long	17" - 20"	111473	0-35664-11473-2
Large	15"-18"	Regular	14" - 17"	111474	0-35664-11474-9
		Long	17" - 20"	111475	0-35664-11475-6
X-Large	18"-20"	Regular	14" - 17"	111476	0-35664-11476-3
		Long	17" - 20"	111477	0-35664-11477-0

For Knee High, measure calf circumference and length from heel to bend of knee.

THIGH HIGH

Size	Calf Circumference 	Length 		Catalog #	UPC Code
Small	9"-12"	Short	25" - 29"	111478	0-35664-11478-7
		Regular	29" - 33"	111479	0-35664-11479-4
		Long	33" - 38"	111480	0-35664-11480-0
Medium	12"-15"	Short	25" - 29"	111481	0-35664-11481-7
		Regular	29" - 33"	111482	0-35664-11482-4
		Long	33" - 38"	111483	0-35664-11483-1
Large	15"-18"	Short	25" - 29"	111484	0-35664-11484-8
		Regular	29" - 33"	111485	0-35664-11485-5
		Long	33" - 38"	111486	0-35664-11486-2
X-Large	18"-20"	Regular	29" - 33"	111487	0-35664-11487-9
		Long	33" - 38"	111488	0-35664-11488-6
XX-Large	18"-21"	Regular	25" - 32"	111489	0-35664-11489-3
		Long	32" - 38"	111490	0-35664-11490-9

For Thigh High, measure calf circumference and length from heel to gluteal fold.



Care instructions:

- Machine wash in lukewarm water with delicate fabric detergent.
- Tumble dry on low heat.
- Do not bleach, dry clean or iron.

Fiber contents:

- Knee High: 85% Nylon; 15% Spandex
- Thigh High: Small – X-Large: 79% Nylon; 17% Spandex, 4% Silicone
XX-Large: 77% Nylon, 19% Spandex, 4% Silicone

NOTE: Please refer to the product label and/or package insert for full instructions on safe use of these products.



BSN medical Inc. 5825 Carnegie Blvd. Charlotte, NC 28209-4633
Tel. (+1) 704 554 9933 Fax (+1) 800 835 4325 www.jobst.com
To order toll-free: JOBST (+1) 800 537 1063

57010 RN

© 2007 BSN medical Inc.

Printed in U.S.A.



Comfort, Health and Style!

REV 03/07